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

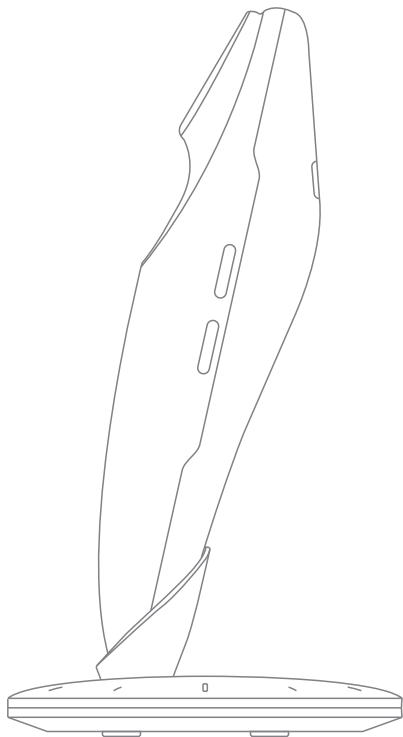
/ Catalogue

01-Main attachments	01
02-Product structure	02
03-Mode and level Display	03
04-Bind APP to use device	04
05-Use Device without binding APP	05
06-About charging	08
07-Product specification	09
08-Safety Caution	09


01 Main attachments

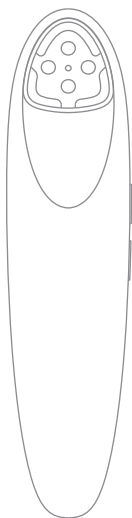
Device

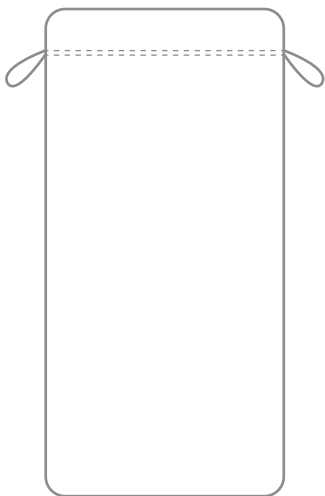Storage Bag


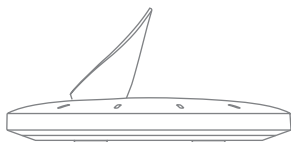
Charging base

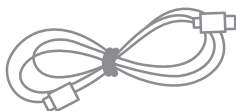USB cable




02 Product structure

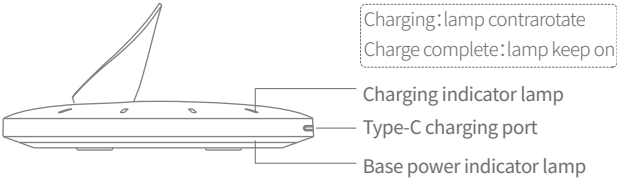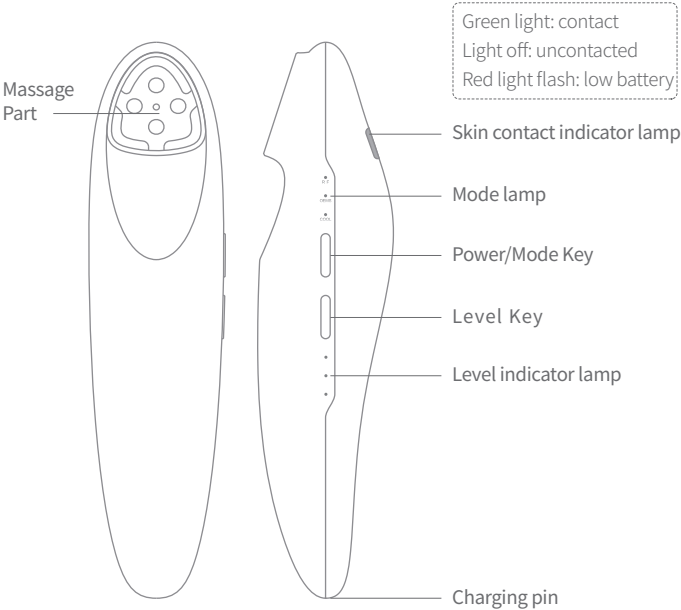Tips:

- ① When using COOL mode, the skin contact indicator is off;
- ② Low power, skin contact indicator lamp flashing red, and the device continuous "beep".

03 Mode and level Display

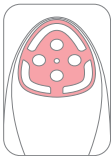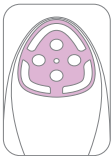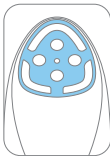Mode	anti-wrinkle (RF)	Tighten Rejuvenation (OEMS)	Soothe and relieve (COOL)
Light display			
Technology application	RF+ Red light	RF+OEMS+ Red blue light	Blue light +Cool
Mode effect	Daily care Anti-winkle & tighten	antiedema, V face Stay up late First Aid	after-sun repair Shrink pores
Single care time	6 min	6 min	3 min
Frequency of use	2-3 times a week		Unlimited of use
Level selected	Level 1、2、3	Level 1、2、3	None
Important notice	<ul style="list-style-type: none"> * Keep face clean; * Cover massage part with Preparation Gel before use; * After using this mode, if you need to use other modes, please wait at least 10 minutes; * After use, clean the massage part with dry cloth or tissue, do not rinse with water; * Do not use with alcohol-based cosmetics. 		<ul style="list-style-type: none"> * Keep face clean;

04 Bind APP to use device


-- Download APP

scan the QR code below or the QR code on the box to download.

-- Device connect APP

- 1. Open the APP and ensure that the Bluetooth of the mobile phone is on;
- 2. According to the instructions in the APP, connect the RF device , enter the care page;
- 3. Start beauty care according to APP guidance.

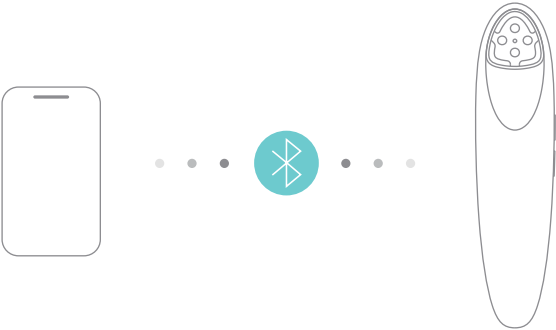
05 Use Device without binding APP


A.RF mode(anti-wrinkle) 、 OEMS mode(tighten) can be used as follows

01 - Cleanse and dry your face.

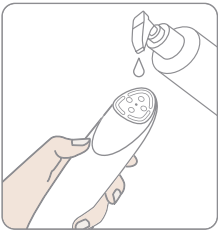02 - Apply appropriate amount of Preparation Gel to the massage p-
art or facial area to be treated. (Pl-
ease add gel appropriately accord-
ing to face sensation during use)

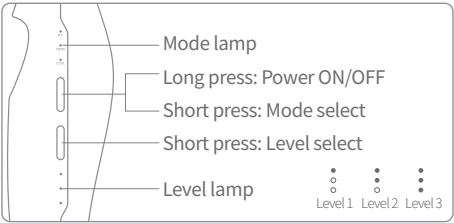
Cation: Be sure to use level 1 for the first time or before getting familiar with it.


03 - Long press the Power/Mode button for 2 seconds to start the device;

Tip: Boot default RF mode, the default level is "level 1", according to your own needs to adjust the mode and level, the indicator will change to the change of mode and level (the first use please adjust level 1 to adapt).

04 - When using, stick the massage part completely on the skin and slowly make a circular motion, as shown in the figure below. The left and right face both for 3 minutes, the device will beep to remind you to switch.

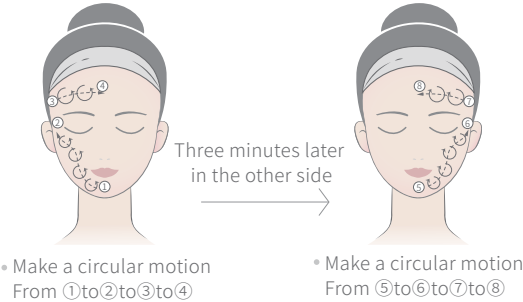Tip: To avoid high temperature scald, do not stick to the same position of the skin for a long time; Do not use on the lips, eyelids, intraorbital bones, thyroid, and cervical lymph nodes.

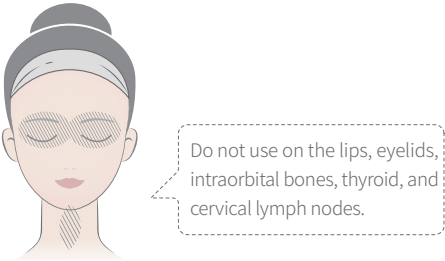05 - 6 minutes to the end of the care, the device automatically power off after beep, the indicator goes out (in the process if you want to stop, long press the power/mode key to power off).

06 - After use

Skin care--The Preparation Gel is not leave-in type, it is recommended to clean the face with water after use, and then begin daily maintenance.

Instrument maintenance--After use, clean the massage part with a dry cloth or tissue, do not rinse with water.


B.COOL mode can be used as follows:

01 - Cleanse and dry your face.

02 - Long press the Power/Mode button for 2 seconds to power on, then short press the button to switch the mode to COOL mode

Tip: Cannot switch levels in COOL mode.

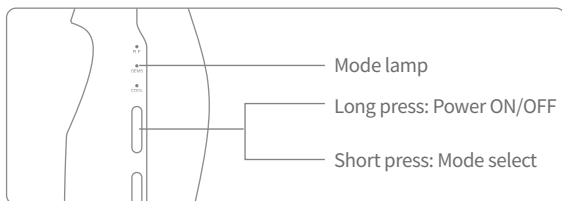03 - During using, apply the massage part to the skin completely, the operation method is shown below. The left and right face press for both 1.5 minutes, the device will beep to remind you to switch;

Apply the massage part to skin completely, Press point one by one.


- Gently massage the designated areas on the face for targeted skincare,
- Left face for about 1.5 min,
- Right face for about 1.5 min,
- Power off after 3 min.

04 - The device will beep and automatically power off after 3 minutes, and the indicator will go off (if you want to stop during use, long press the power/mode key to shutdown).

06 About charging

* Charging method

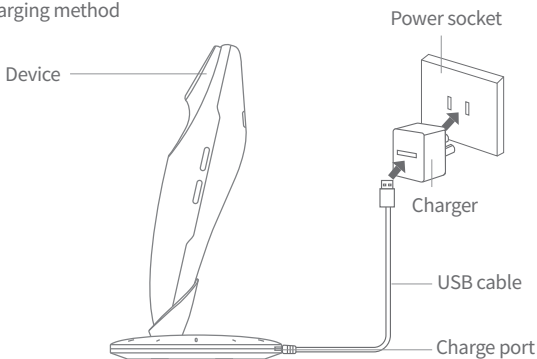* There may be differences between the picture and the real thing. Please refer to the real thing For power saving and safety, please unplug the charger and the USB cable after charging.

* Charging state

No charge: the base power indicator is steady on, and the charging indicator is off;
Charging: the base power indicator is steady on, and the charging indicator turns counterclockwise;

Fully charged: The base power indicator is steady on, and the charging indicator is all on.

* Low power display

The skin contact indicator lamp flashes red, and the device continuous beep.

* The instrument can be turned on during charging.

07 Product specification

Name	Skin Care Device (OGP Smart Skin Renewal RF Device)
Type	OGP collagen 3
Input	DC5V \Rightarrow 2.0A
Power Consumption	\approx 16W
Max RF Output	200V \sim 1MHz
Charge Time	about 3 hours *Charge time depend on the charging environment.
Electrode	Φ =5mm*4
Standards	IEC60335-1,IEC60335-2-23, IEC60335-2-115
Usage Notice	Please leave more than 24 hours between treatments
Place Of Origin	China

08 Safety Caution

Please read this safety cautions carefully before using. The cautions are to ensure the correct use of the product and to prevent personal and property losses of users and others. Please be sure to follow them.

> Signs description

The following signs are used to classify and describe the extent of danger, injury, and property damage resulting from disregarding sign instructions and improper use.

DANGER

Indicates a danger that is likely to result in serious injury or death.

WARNING

Indicates a potential danger that could result in serious injury or death.

CAUTION

Indicates that damage may be caused to the home, property, poultry, or pets.

The following signs are used to classify and describe of instructions that need to be followed.

 The sign is used to alert the user to prohibited operations.

 The sign is used to remind users of the steps they must follow to operate the product safely.

 **DANGER- Prevent accidents**

 Users of the following electronic medical devices and home therapy devices are not allowed to use this product.

- Implantable medical devices such as cardiac pacemakers
 - Life-sustaining medical instruments such as artificial heart and lung
 - Wearable medical instruments such as electrocardiograph and drip controller
 - Failure to do so may result in misoperation of electronic medical devices and home treatment devices, which may cause accidents or physical discomfort.
-

 **WARNING - Prevent accidents**

 To prevent skin and physical discomfort, please do not use this product:

- Heart patients
 - People with facial nerve injury
 - Diabetic patients
 - Patients with hemorrhagic disease, hemophilia and other blood abnormalities
 - People with hypertension
 - People undergoing dental treatment
 - People with exfoliated retina or who have had exfoliated retina in the past
 - People with high myopia (above -6D)
 - People who have received metal, silicone, plastic implants on their face due to fracture, plastic, etc
 - Young and growing children, people unable to express their wishes clearly
 - People with sequelae of head and eye due to traffic accidents, etc
 - People who have undergone eye surgery (for those who have undergone LASIK surgery, please consult your primary physician)
 - Otherwise it could cause an accident or physical discomfort.
-

WARNING

The following persons should not use this product:

- People with skin diseases and allergic dermatitis
- Those who are pregnant, menstruating or have just given birth
- People with slow heat sensitivity and insensitive skin are forbidden to use this product
- Otherwise, it may cause skin problems, physical discomfort or accidents.

Do not use this product in the following situations:

- When feeling unwell
- When you feel particularly tired during drinking or taking medicine
- When there are cuts, eczema and lumps on the skin
- When you have skin diseases such as dermatitis and sunburn
- Otherwise it may cause skin or body discomfort.

⊗ Do not use this product in the following situations:

- The site was depilated on the same day
- The site has received plastic surgery
- Otherwise it may cause skin or body discomfort.

⊗ Do not use this product in the following areas:

- Lips, eyelids, intraorbital bone, thyroid gland, and cervical lymph nodes
- Otherwise it may cause accidents or physical discomfort.

! Please consult your doctor before using this product for the following groups:

- People with sensitive skin
- Long-term severe patients
- Patients with acute diseases
- Patients with febrile diseases
- Patients with infectious diseases
- Patients with malignant tumors
- Those who have received treatment that requires a certain amount of recovery time at beauty salons, clinics, etc., or have received treatment that causes severe skin redness and swelling within the last 2 weeks
- Otherwise, it may cause an accident or cause skin or physical discomfort.

! Please stop using this product immediately if you feel any pain or abnormality on your skin or body.

- Otherwise, it may cause an accident or cause skin or physical discomfort.

! Do not allow children to use it or avoid disabled people to use it alone. Minors should obtain the consent of their guardians and use it with them.

- Failure to do so could cause an accident or injury.

! Do not allow use by persons who are unable to express themselves and have very insensitive skin perception.

- Failure to do so could cause an accident or injury.

! Keep the host /AC adapter out of reach of children. Do not allow children to use or play with it.

- Failure to do so could cause an accident or injury.

! Please properly handle the packaging of this product and keep it away from infants and children.

- Otherwise, accidents such as suffocation may occur.

> The product

⊘ Do not use this product near flammable materials.

· Alcohol, gasoline, thinner, spray, nail polish remover, etc.

- Failure to do so could result in fire or injury.

⊘ Do not put into the fire or heat. Do not charge, use or place this product under the following circumstances:

· Places with heavy moisture such as bathrooms; Dusty places; A place in direct sunlight; Places prone to water; Places prone to tilting, shaking or being hit; Outdoors; A place with strong electrical or magnetic waves.

· Places accessible to children; Places with soot and steam; High temperature place; Around heating appliances (stoves, etc.); A place with a large temperature difference (next to air conditioning, heating, etc.).

- Failure to do so may result in fire, electric shock, accidents, breakdowns and falls.

⊘ **Do not modify or repair the product yourself.**

- Failure to do so could result in fire, electric shock, or injury.
- If you need maintenance, please contact customer service (such as battery replacement).

⊘ **Do not disassemble the product except discard it.**

- Failure to do so may result in electric shock, injury, fire or instrument damage.

⚠ **WARNING**

⊘ **Keep this product out of the reach of children and infants. In addition, do not allow them to touch or handle this product.**

- Otherwise accidents and injuries could result.

> Note the following items

⊘ **Please prevent the product from falling or being impacted.**

- Otherwise, the instrument may be damaged.

⊘ **Do not use this product in the same area for a long time.**

- Otherwise, it can cause skin damage.

❗ **Before use, please check the probe surface for cracks, notches, deformation and other damage.**

- Otherwise, skin injury may occur.

> Power supply

⊘ **Do not connect the adapter to or remove the adapter from a home electrical outlet with wet hands.**

- Electrocuting or injury may occur.

- ⊘ Do not immerse the body in water or clean it with water.
 - This could result in electric shock or fire.
- ⊘ Do not use or store in places with high humidity, such as the bathroom. Do not store in places where prone to water accumulation (such as the sink).
 - Electrocutation or fire may result.
- ⊘ Do not exceed the rating of household electrical outlets or connections when using.
 - If too many plugs are attached to an electrical outlet and it exceeds its rating, it may cause a fire due to overheating.
- ❗ Be sure to use 5.0 v=2.0 A specifications of the power adapter to charge to our products.
 - Failure to charge or damage to the instrument may result.
- ❗ Be sure to use 100-240V.
 - Failure to do so could result in fire or electric shock.
- ❗ Please fully insert the adapter or power plug.
 - Failure to do so may result in fire or electric shock.
- ⚠ WARNING
 - ❗ For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance.
 - ❗ Be sure to remove the adapter from your home electrical outlet while cleaning.
 - Otherwise it may result in electric shock or injury.
 - ❗ When the main engine and charging end are wet, please thoroughly dry the water.
 - Electric shock, short circuit and fire if charged with residual moisture.
- ⚠ CAUTION
 - ⊘ Do not allow the needle or foreign object to attach to the power plug and power cord plug.
 - Otherwise, it may cause electric shock or fire due to short circuit.

⊗ Do not wrap the power cord tightly around the adapter during storage.

- Otherwise, the power cord may be broken and an electric shock may occur or a fire may be caused by a short circuit.

❗ Disconnect the power adapter from your home electrical outlet when not charging. When removing the power adapter, hold the adapter, not only the power cable.

- Otherwise, electrical leakage may result from the deterioration of insulation performance, which may lead to electric shock or fire.

> If an exception or failure occurs

❗ If abnormal or malfunction occurs, please stop using and unplug the power plug immediately.

- Otherwise it may result in electric shock, injury or fire.

❗ Abnormal or malfunctioning condition:


- Deformation or abnormal heating of the body, adapter or power cord;
- The body, adapter or power cord emit a burning odor;
- Abnormal sound from the body, adapter or power cord during use or charging.

> The abandonment of the built-in rechargeable battery


- this device use rechargeable lithium battery, please according to local regulations for scrap.
- before discarded appliances, batteries must be removed from the appliance.
- take out the battery, the instruments must be without electricity.
- removing the battery, you need to use the lever and a screwdriver, please pay attention to safety to avoid injury.
- removing the battery, please avoid hard objects to cut the battery. * May cause battery heat, fire, etc.

> Operation


Remove the battery by referring to steps ①~④ in the following figure.

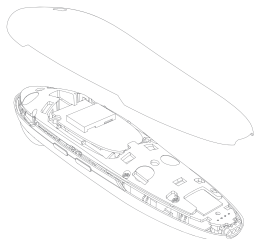
① Use the pry bar tool to pry open the rear shell

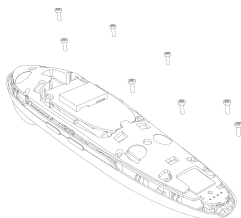② remove the screws of the rear shell support

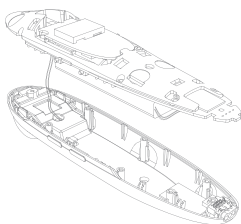③ remove the support

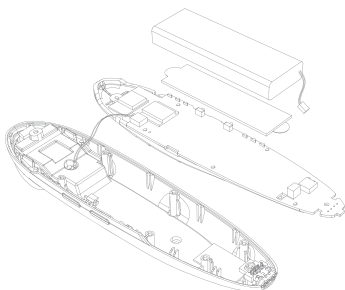④ remove the power cable, remove the instant glue, and remove the battery

⚠ CAUTION

- use this guide applies only to scrap the host, and built-in rechargeable lithium batteries need operation guidelines in accordance with the relevant laws and regulations on recycling.
- this user guide can't be used in normal product maintenance.
- if removed resulting in failure of the host, please will can't provide a normal after-sales service.

> Remove the battery when the product is discarded

DANGER

This rechargeable battery is only used for this product. Do not use it for other products.

- Do not charge the battery that has been removed from the product.
- Do not heat it or put it in the fire.
- Do not knock, disassemble, modify, or scratch with your nails.

Do not connect the positive and negative poles of the battery with metal objects.

- Do not carry or store batteries with metal accessories such as necklaces and hair clips.
- Do not use or place batteries in high temperature environments (such as direct sunlight or near other heat sources).
- Do not tear off the battery skin.
- Otherwise, it will cause overheating, fire or explosion.

WARNING

Once the rechargeable battery has been removed, do not place it within the reach of children and infants.

- Accidental ingestion of batteries can result in personal injury, if this occurs, seek medical attention immediately.

If the electrolyte leaks and contacts the eyes, please do not rub the eyes, please rinse thoroughly with water immediately.

- This may lead to eye injury, if this occurs, seek medical attention immediately.

CAUTION

If the electrolyte leaks and contact skin or clothing, rinse immediately with water.

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Manufacturer: HangZhou SHIGUANGJI Intelligent Electronic Technology Co.,Ltd

Address: Room 617-620, Building 5, No.2 kejiyuan Road, Baiyang Street, Qiantang District, Hangzhou City, Zhejiang Province.